

Supermini200

SEQUENTIAL BENCHTOP WDXRF SPECTROMETER

ELEMENTAL ANALYSIS BY
WAVELENGTH DISPERSIVE
X-RAY FLUORESCENCE



Rigaku

POWERING NEW PERSPECTIVES



TABLE OF CONTENTS

At About ½ the Price and ¼ the Size, the Supermini200 Changes Everything	3
It Is Time for You to Reconsider WDXRF.....	4
How WDXRF Works.....	4
Newly Designed and Simplified Software	4
High-Performance Elemental Analysis	5
Rich Hardware Features Ensure Sample Measurement Flexibility.....	6
Specifications	7

AT ABOUT ½ THE PRICE AND ¼ THE SIZE, THE Supermini200 CHANGES EVERYTHING

Wavelength dispersive X-ray fluorescence (WDXRF) is one of the most powerful and well-established techniques for elemental analysis. It has several advantages, including light element sensitivity, exceptional elemental resolving power, and low limits of detection. Historical drawbacks of the WDXRF technique have been that such instruments are typically large, expensive, and require special utilities. With a size that is a quarter of traditional WDXRF units and a price that is approximately half of floor standing units, the Supermini200 changes everything.


The X-ray fluorescence (XRF) experiment begins with exposing a sample to high-energy photons from an X-ray tube, which induces transitions of electrons between atomic orbitals and results in the emission of fluorescent photons. By measuring the energy and

intensity (count rate) of these photons, qualitative and quantitative information about the elemental composition is obtained.

Measurement quality, usually judged as precision (repeatability), is a function of many factors. Two of the primary metrics are elemental peak resolution and photon counting rate.

- **Higher resolution results in narrower and taller peaks, which improves precision and sensitivity for elements with similar fluorescent peak energies**
- **Higher counting rates improve precision and/or shorten the overall analysis time**

Analyzing samples that contain only a few non-overlapping heavy elements is easy for any XRF instrument. Measuring light elements is always more difficult since the fluorescent yield of emitted photons is proportional to the fourth power of the atomic number (Z^4). Lighter elements not only emit fluorescent photons with lower efficiency, but their peaks can also overlap with each other as well as with the L and M emission lines of heavier elements. The Supermini200 is a WDXRF spectrometer and has much better spectral resolution than EDXRF systems. Therefore, line overlaps in the Supermini200 are minimized.

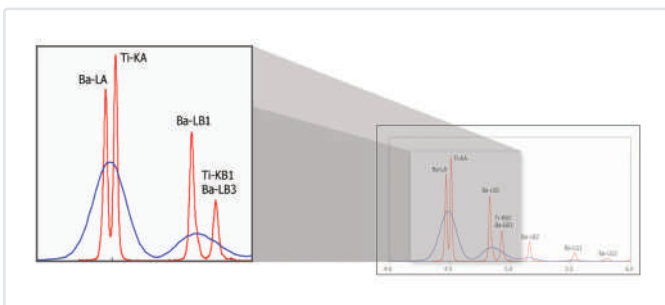

MINING AND ORES

Suitable for the analysis of various kinds of ores and concentrates for base and rare metals, such as copper, nickel, iron, and coal. The Supermini200 is ideal for use in small laboratories located at mining sites, satellite laboratories, or central laboratories as a backup system for large WDXRF.

CEMENT AND RAW MEAL

High precision, sensitivity for light elements, and compliance with ASTM C114 make the Supermini200 an ideal tool for the cement industry as a primary or backup analyzer for raw materials, raw meal, clinker, and finished cement.

METALS AND ALLOYS

The WDXRF method of the Supermini200 is the preferred technique for analyzing steel-making slags, steels, specialty alloys, aerospace materials, and turbine blades as a total process control tool.



LUBRICANTS AND ADDITIVES

Analysis of lubricants and additives for Ca, Cl, Mg, P, S, Zn and Mo is no challenge for the Supermini200, in spite of peak overlaps and the presence of light elements.



IT IS TIME FOR YOU TO RECONSIDER WDXRF

Rigaku's Supermini200 combines all the advantages of traditional WDXRF elemental analysis systems into a smaller, more economical package. Now is the time to reconsider WDXRF as one of your primary elemental analysis techniques. Compared to other techniques, such as ICP and AA, there is no chemical preparation step required, making sample preparation and cleanup much simpler.

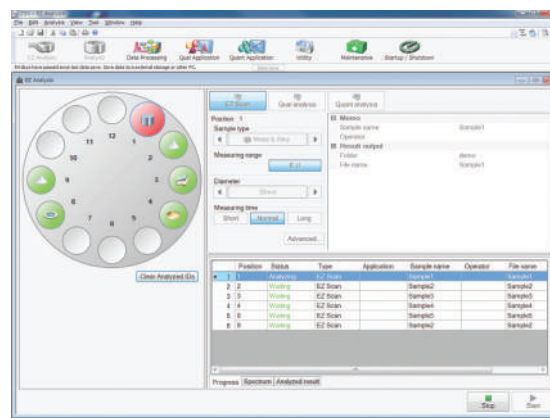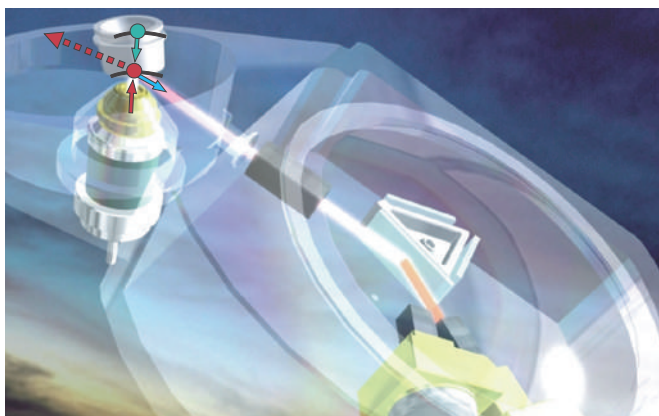

HOW WDXRF WORKS

Excitation

An X-ray photon from the tube removes one of the inner-shell electrons and excites the atom.

Transition

Another electron from a higher-energy orbital transitions down to fill the lower "vacant" orbital.

Fluorescence

The difference between orbital energies is emitted as an X-ray photon characteristic to the element.

Dispersion

X-rays are physically separated—or dispersed—using crystals that are designed to reflect a different wavelength at each incident angle.

Counting

At each angle, the reflected beam enters a high-speed detector that can count >1M photons/second, providing the capability of achieving very high precision.


NEWLY DESIGNED AND SIMPLIFIED SOFTWARE

EZ Analysis is a new feature that makes everyday routine operation simpler. A single interface contains everything you need to know about the status of your samples, what your data measurement parameters are, what your data analysis parameters are, and a running output of your results.

EZ Scan allows you to analyze unknown samples without any prior setup and with only a few clicks of the mouse to get started. Combined with Rigaku's SQX fundamental parameters software, it provides the most accurate and rapid XRF results possible. SQX is capable of automatically correcting for all matrix effects, including line overlaps. SQX can also correct for various atmospheres, impurities and sample film absorption. Increased accuracy is achieved using a matching library. At each angle, the reflected beam enters a high-speed detector that can count >1M photons/second, providing the capability of achieving very high precision.



Air-cooled X-ray Tube

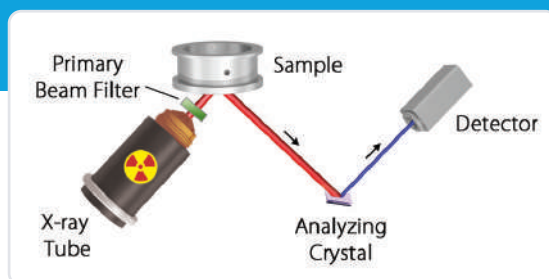

HIGH-PERFORMANCE ELEMENTAL ANALYSIS

UNIQUE ADVANTAGE #1

The Supermini200 is the only instrument to give the light element sensitivity and energy resolution of a WDXRF spectrometer with the footprint of an EDXRF spectrometer.

UNIQUE ADVANTAGE #2

A small 200 W air-cooled X-ray tube (pictured above) provides an excellent balance between power, instrument size requirements, and utility usage.

UNIQUE ADVANTAGE #3

The primary beam filter improves specific elemental analyses: Zr filter for Cd, Ru, Rh, Pd, Ag, or In analysis or optionally an Al filter for trace K analysis and for Ag on film.

UNIQUE ADVANTAGE #4

The Supermini200 employs a scintillation counter (SC) for heavy elements and a gas-flow proportional counter (F-PC) for light elements. This ensures high-sensitivity and high-speed counting.

UNIQUE ADVANTAGE #5

The Supermini200 software has been rewritten to emphasize ease-of-use, making it especially easy to learn for new operators and improving efficiency for experienced users.

UNIQUE ADVANTAGE #6

Samples can be added or removed from the 12-position automatic sample changer without halting the current data measurement.

RICH HARDWARE FEATURES ENSURE SAMPLE MEASUREMENT FLEXIBILITY

Solid samples, such as metals, alloys, and pressed powder briquettes, are measured under vacuum. Liquid samples, such as aqueous solutions or oil, are poured into liquid cells, covered with sample film and measured under helium. A 12-sample changer is standard, enabling operators to carry out routine analysis smoothly. For added operator flexibility, during measurement you can replace samples on the changer without interrupting the measurement.



Solid Sample Cup



Automatic Sample Changer



Liquid Sample Cup

The Supermini200 is equipped with a three-crystal exchanger with LiF(200) and PET mounted as standard crystals. RX26, RX9, or Ge can be added optionally.

Crystal		Atomic Number											
		1	10	20	30	40	50	60	70	80	90		
Standard	LiF(200)			²² Ti									⁹² U
	PET		¹³ Al		²¹ Sc								
Optional	RX26		⁸ O		¹² Mg								
	RX9			¹⁵ P		¹⁷ Cl							
	Ge			¹⁵ P		²¹ Sc							

Crystal Selection and Functionality

PETROLEUM AND BIOFUELS

Ultra-low detection limits for S, P, and Cl, together with the ability to analyze many other elements, make the Supermini200 an attractive analytical tool for petroleum refineries and biofuel plants.

LAYERS AND THIN FILMS

Researchers and process engineers utilize the Supermini200 for rapid, non-destructive analysis of composition and thickness of thin films, including multilayer structures such as photovoltaic cells, using Rigaku's advanced thin film FP software.

COATINGS

The Supermini200 can determine coating weight and elemental composition of the metallic and non-metallic coating layer, as well as surface treatment on steel or aluminum sheets.

ENVIRONMENT


Plastics with toxic additives, incinerator sludges, and contaminated soils can be analyzed thanks to Rigaku's powerful SQX FP semi-quantitative software.

SPECIFICATIONS

Specifications	
Elemental Coverage	Oxygen (O) through uranium (U)
X-ray Tube	Pd target
X-ray Power	200 W
Generator	50 kV, 4 mA
Cooling	Air cooled
Beam Filter	Programmable Zr standard Al optional
Crystals	LiF(200) and PET standard RX26, RX9, or Ge optional
Detectors	Flow proportional counter Scintillation counter
Sample Size	10-position turret: 52 mm (2.0 in) diameter maximum 33 mm (1.3 in) height maximum 12-position turret: 44 mm (1.7 in) diameter maximum 33 mm (1.3 in) height maximum
Auto-sampler	10- or 12-position turret
Sample Spinner	Standard, 30 rpm
Atmosphere	Vacuum standard Helium optional
Vacuum Pump	Rotary pump
Helium Flush	15 – 65 PSIG, 0.5 l/min
Environment	15 – 28°C (59 – 82°F) temperature <75% relative humidity
Power Supply	100 – 120V (50/60 Hz) 15A or 200 – 240V (50/60 Hz) 10A

Software	
Operating System	Windows® 11
Options	SQX software with FP Matching library SQX scatter FP method Fused bead correction Line Overlap Correction using Theoretical Intensities (LOCTI) Quant scatter FP

Dimensions and Mass	
Width	580 mm (22.8 in)
Depth	680 mm (26.8 in)
Height	670 mm (26.4 in)
Weight	100 kg (220.5 lb)
Vacuum Pump	170 mm (W) x 500 mm (D) x 310 mm (H), 28 kg (6.7 in x 19.7 in x 12.2 in, 61.7 lb)

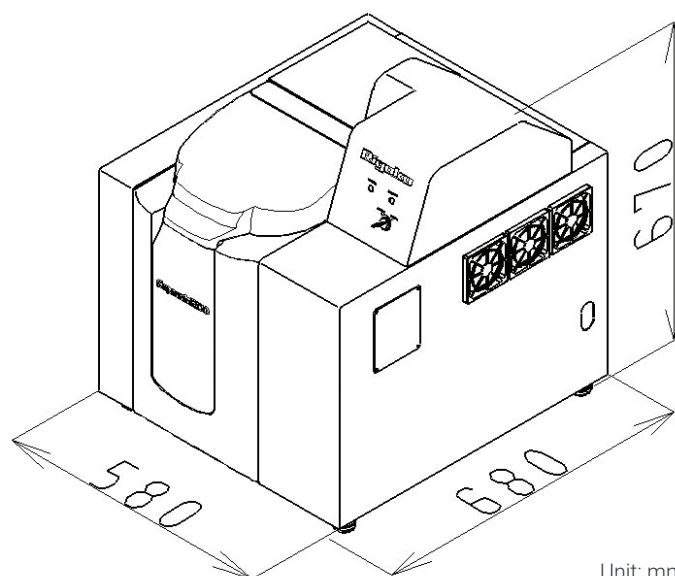

Unit: mm
Weight: 100 kg

Supermini200

SEQUENTIAL BENCHTOP WDXRF SPECTROMETER

Rigaku Corporation and its Global Subsidiaries

info@rigaku.com | www.rigaku.com

rigaku.com/products/wdxrf



ELEMENTAL ANALYSIS BY WAVELENGTH
DISPERSIVE X-RAY FLUORESCENCE

AXEL

NITTO KAGAKU CO., Ltd.

1-300-02 Automatic Mortar ANM-200



Features

- It has powerful pulverization capability by automatic rotation of a pestle and a mortar.
- It doesn't take a place by compact design but carrying is easy.
- It is an omnipotent type by adapter exchange, such as crockery, alumina and agate.

Spec

- Size (mm): 210 x 350 x 600 - 775
- Motor: Upper and lower 1 unit each 100VAC 46W/41W
- Power cord length: 4.5m (with 2P plug)
- Rotational speed (50/60Hz): Pestle unit/100/120rpm, Mortar/6/7rpm
- Weight: 18kg
- Model number: ANM-200
- Accessories: Mortar AN-20, Pestle AL-20B, Spatula AN-20H

***This product might require extra fees.**

▶ Package size : 270×420×690 mm 20 kg

AXEL

AS ONE Corporation

6-522-01 Spoon (Stainless Steel) 3 Set (150, 165, 180mm)



匙

3本組 (150mm・165mm・180mm)

Spec

- Material: Stainless steel (SUS410)
- Three-piece set (150mm, 165mm, 180mm)
- Type: Spoon
- RoHS3 (EU)2015/863 Certificates can be issued

▶ Package size : 40×200×10 mm 60 g

Manual Hydraulic Press - for FTIR & XRF Sample Preparation

The 15 ton and 25 ton Manual Hydraulic Presses have been designed to handle a wide variety of pressing applications. They are especially suited for the preparation of solid samples for XRF analysis and KBr discs for FTIR analysis Specac Evacuatable Pellet Dies. The Presses can also be used with Specac Heated Platens for applications such as the preparation of thin polymer films.

Key Features of the Manual Hydraulic Press

- Rugged hardwearing design
- Polycarbonate safety guards
- Adjustable upper bolster
- Adjustable pressure relief valve
- Vacuum Port
- Low Pressure Gauge option



Ordering Information

Manual Hydraulic Press

GS15011	15T Manual Hydraulic Press
GS25011	25T Manual Hydraulic Press

Spares and Consumables

GS15100	Replacement Seals and Gaskets kit for Manual Hydraulic Presses
GS15101	Hydraulic Oil (1litre)

Spares and Consumables

GS15051	0-1 Ton Gauge Conversion Kit
GS15052	0-2 Ton Gauge Conversion Kit
GS15055	0-5 Ton Gauge Conversion Kit

Specifications

Max height (at handle):	610 mm
Width:	310 mm
Depth:	190 mm
Weight:	50 kg
Lower piston stroke:	25.4 mm
Upper bolster screw travel:	89 mm
Distance between pressing faces (max):	152 mm
Distance between pressing faces (min):	38 mm
Lower pressing face diameter:	86 mm
Upper pressing face diameter:	32 mm
Max width of sampling area (side to side):	143 mm
Max depth of sampling area (front to back):	141 mm

United Kingdom
sales@specac.co.uk
+44 (0) 1689 892 902

United States
sales@specac.com
+1 866 726 1126

China
frank.li@specac.com

Singapore
kamhar.woo@specac.com



AXEL

3-7534-11 KBr Plate KBr Cutting STJ-0204S



Spec

- Piece Size (mm): 3 x 3 x 0.75
- Quantity : 1 piece (50 pieces)
- material: KBr

▶ Package size : 35×70×30 mm 40 g

AS ONE offers everything you need from laboratory instruments and consumables, factory MRO, to hospital & nursing care supplies



AXEL

Potassium Bromide Premium Quality 25g



Spec:

Chemical Formula: KBr

Standard: Premium Quality (JIS Certified Reagents which meet the premium standard)


Concentration and Purity: 99%

Molecular Weight: 119

Sales Unit: 25g/bottle

AXEL

4-3727-03-30 Agate mortar (deep) ϕ 88.9mm POMP 306

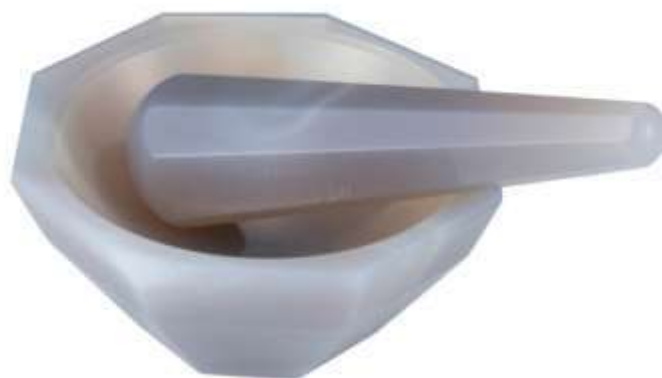

Spec

- Outer diameter (mm): ϕ 88.9
- Number of pieces: 1
- Material: Agate
- Accessories: Pestle

AXEL

AS ONE Corporation

6-522-01 Spoon (Stainless Steel) 3 Set (150, 165, 180mm)



匙

3本組 (150mm・165mm・180mm)

Spec

- Material: Stainless steel (SUS410)
- Three-piece set (150mm, 165mm, 180mm)
- Type: Spoon
- RoHS3 (EU)2015/863 Certificates can be issued

▶ Package size : 40×200×10 mm 60 g

AXEL

AS ONE Corporation

1-2941-02 Hot Plate HHP-170D



Features

- Possible to accurate temperature control by PID control method.

Spec

- Model number: HHP-170D
- Highest temperature: 300 °C
- Plate size: 170 x 170mm
- Plate material: Aluminum
- Digital
- Temperature control system: PID control
- Temperature settings unit: -10°C
- Heater capacity: 600W
- Material: Body/stainless steel (SUS304), steel (SPCC)
- Safety device: Over current fuse, internal temperature fuse, high temperature warning display function
- Power supply: 100VAC 50/60Hz
- Main unit size: 170 x 200 x 125mm
- Weight: 3kg
- RoHS3 (EU)2015/863 Certificates can be issued

► Package size : 300×235×170 mm 3.3 kg

AXEL

AS ONE Corporation

1-001-01 Desiccator 300 x 345 x 525mm LH



Features

- By form of lid is uniform and it has very high airtightness.
- Shelf receptacle is arbitrarily movable and enters easily from small to big thing.


Spec

- Material: Body/transparent PMMA (acrylic), Shelf board/PS (polystyrene), Rubber foot/natural rubber
- Frontage size: 258 x 458mm
- Shelf board dimensions: 257 x 270 x 4mm
- Shelf board load capacity: Approximately 2kg/board
- Total withstand load: 10kg
- Shelf board mounting pitch: 1 stage in the bottom/50mm, Other/30mm (minimum shelf board spacing 20mm)
- Weight (including shelf board): 4.5kg
- Accessories: Shelf board x 3, silica gel plate x 1
- Outer dimensions: 300 x 345 x 525mm
- Inner dimension: 285 x 275 x 485mm
- Model number: LH

► Package size : 390×380×635 mm 6.2 kg



1-3122-01 Dehumidifying Agent 10 Pcs M-50FR

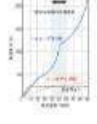

Features

- Max.moisture absorption volume is 200% of weight of itself.
- Effect lasts for long time also at low humidity.

Spec

- Weight: 50 g / unit
- Model number: M-50 FR
- Size (mm): 115 x 138
- Quantity: 1 case (10 pieces)
- * It cannot be played back.
- RoHS3 (EU)2015/863 Certificates can be issued

▶ Package size : 290×290×55 mm 600 g